

Features

1 Three-way syringe (air/water/spray) with sterilisable tips (hose and tips detachable)

1 Electric micromotor with LED light and 30 programs (0–40,000)

1 Electric scaler with LED light

Instrument hoses detachable

Spray adjustment regulators for each instrument

Integrated spray water system for syringe and instruments (check valve prevents the back flow of water)

1 Spray water tank: 0.75 litres

1 High-performance suction hose 300 l/min. Continuously adjustable. Unique BPR Swiss technology
Adapter for saliva ejector tips and oral-dry-cup

1 Suction fluid container: 1 litre (automatic overflow protection)

Pressure regulator-filter

Sturdy, air- and water-tight case, will not scratch, dent or corrode

Detachable front cover

Disc type foot control

Integrated silent compressor – oil-free: Max. Pressure: 4 bar

DENTA-TROLLEY with drawers for instruments and consumables

Technical Data

Operating pressure: 3–4 bar

Each instrument is individually adjustable

Weight SMART-PORT Premium: 12 kg

Weight DENTA-TROLLEY: 8 kg (space for instruments and consumables)

Dimensions SMART-PORT Premium: 380mm × 205mm × 530mm

Dimensions DENTA-TROLLEY: 410mm × 400mm × 550mm

Warranty: 2 years

Options: Contra-angles, air polisher module, additional electrical motor, DENTA-CASE etc.

Smart-Port Premium Comfort made easy

The cleverly designed SMART-PORT Premium is the perfect device for the mobile dentist.


In an 12-kilogram light case, SMART-PORT Premium combines the most important functions and operations of a dental treatment unit and leaves nothing to be desired: Convenient transport, easy setup and dismantling and full performance in continuous operation.

The integrated electric micromotor with LED light features a control unit with 30 selectable programs. Speed, torque, and the selected contra-angle can be freely adjusted and saved.





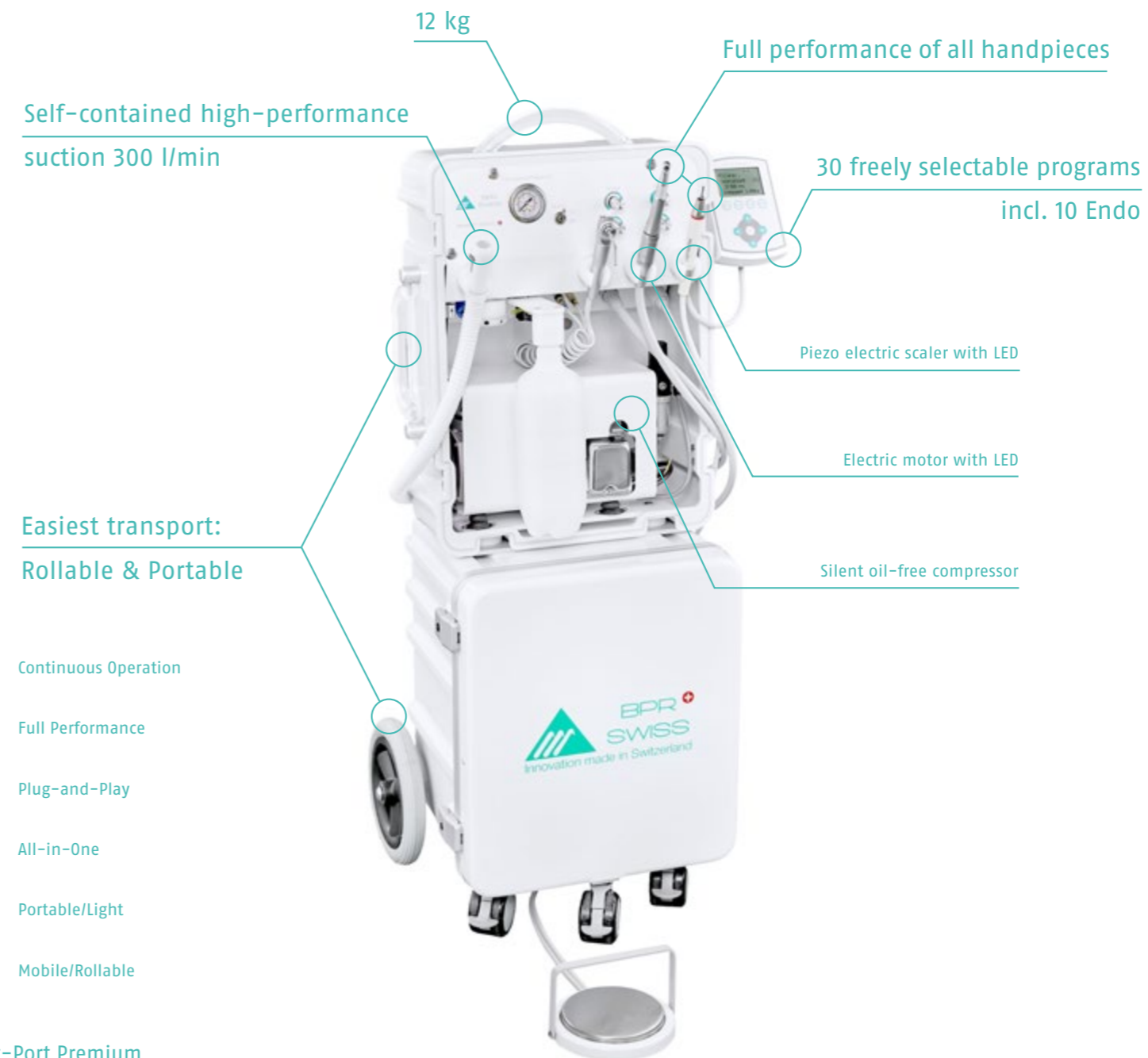

Through the use of various contra-angles, speeds of between 0–200,000 rpm can be achieved. Both the electric motor and the integrated electric scaler with LED light work at their full capacity.

The high performance of SMART-PORT Premium is rounded off by the highly efficient worldwide unique BPR Swiss suction (300 l/min!). The clever design allows for quick and easy assembly and disassembly of the unit. Additionally, the SMART-PORT Premium impresses with its quiet oil-free compressor.

Plug-and-play portable dental units



Easiest transport:
Rollable & Portable

-  Continuous Operation
-  Full Performance
-  Plug-and-Play
-  All-in-One
-  Portable/Light
-  Mobile/Rollable

Smart-Port Premium



4 drawers for instruments/
consumables
Separation red and green zone



On a single trolley transport of
unit/instruments/consumables



Rollable and portable: In case
of obstacles or for stowage the
unit is undocked

Smart-Port Premium